Potenza RETRACTTA

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PHS DO BRASIL LTDA

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Potenza RETRACTTA

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PHS DO BRASIL LTDA

Guilherme Boldt Street, 210, District: Zona Industrial Norte, City: Joinville,

State: Santa Catarina, Country: Brazil - Zip Code: 89239-260

CNPJ: 17.910.389/0001-05 - BRAZILIAN INDUSTRY

E-MAIL: contato@phsdobrasil.com.br

Health Ministry's Operating Licence: X401M5619X7L (8.11992.7)

Notification at ANVISA number: 81199270036

Product Commercial Name: Potenza Retractta

Alternative Commercial Name: Gingival Retraction Cords

Technical Product Name: Restorative dentistry instrument

READ ALL THE INFORMATION IN THIS INSTRUCTION FOR USE CAREFULLY BEFORE APPLYING THE PRODUCT.

KEEP THIS INSTRUCTION UNTIL THERE IS NO FURTHER INTERACTION OF THE PRODUCT WITH YOUR PATIENT.

Technical Support

0800-000-3780

www.phsgroup.com.br

BASIC COMPOSITION:

Cord #000 and #00 – Pima cotton and polyamide.

Cord #0 and #1 – Pima cotton and polyester.

PRODUCT DESCRIPTION:

POTENZA RETRACTTA is a knitted gingival retraction cord developed with a textile structure that combines Pima cotton with polyamide or polyester fibers, depending on the product presentation, combining the absorption properties of natural fibers with the mechanical characteristics of synthetic fibers. This combination provides high mechanical strength, greater structural stability, and a lower tendency to fray during insertion and removal, promoting precise handling and maintaining cord integrity throughout the clinical procedure. The incorporation of synthetic fibers also helps reduce excessive moisture retention within the cord structure, contributing to better control of the operative field during restorative procedures.

Available in different thicknesses identified by color, POTENZA RETRACTTA adapts to different gingival sulcus anatomies and clinical procedure requirements and may be used dry or previously soaked in a hemostatic or astringent solution, according to the technique adopted by the dentist.

PRODUCT' INDICATION:

Potenza Retractta is intended for temporary mechanical retraction of gingival tissue to promote gingival sulcus displacement during restorative and prosthetic dental procedures, such as impressions, crown fabrication, veneers, indirect restorations, and other procedures requiring improved access, visibility, and control of the operative field.

The product may be used alone or in combination with hemostatic/astringent solutions, according to the technique adopted by the professional.

PRECAUTIONS:

- Before use, inspect the integrity of the packaging and ensure that the product is within its expiration date.
- Select the POTENZA RETRACTTA thickness that is most appropriate for the anatomical characteristics of the gingival sulcus and the clinical procedure to be performed.
- Use appropriate instruments for inserting the retraction cord, without sharp edges, in order to avoid damage to the product and gingival tissues.
- Carefully insert the cord into the gingival sulcus, avoiding excessive compression of the tissues and respecting the patient's anatomy.

- The time the cord remains in the gingival sulcus should be limited to the period necessary to achieve the desired retraction, according to the clinical technique adopted.
- When used in combination with hemostatic or astringent solutions, strictly follow the instructions for use and recommendations provided by the manufacturer of the solution used.
- Before removal, when necessary, lightly moisten the cord with water to facilitate its removal and minimize the risk of trauma to the gingival tissue.
- Ensure complete removal of the cord at the end of the procedure, verifying that no fragments or fibers remain in the gingival sulcus.
- After use, dispose of the product in accordance with biosafety standards and the applicable legislation for healthcare waste.

CONTRAINDICATIONS:

POTENZA RETRACTTA is contraindicated in the following situations

- Patients with known hypersensitivity to any of the materials composing the product.
- Presence of acute gingival inflammation, infection, abscess, ulceration, or other soft tissue lesions that may be aggravated by mechanical retraction, unless the procedure is expressly indicated and performed at the clinical discretion of the dentist.
- Patients with uncontrolled periodontal disease when insertion of the retraction cord may worsen the clinical condition.
- Situations in which mechanical gingival tissue retraction is not indicated by the dentist.

POSSIBLE SIDE EFFECTS:

During or after the use of POTENZA RETRACTTA, the following transient effects may occur:

- Gingival discomfort or sensitivity during insertion or removal of the cord;
- Mild gingival bleeding, especially in inflamed tissues or during clinical manipulation;
- Temporary gingival irritation or inflammation;
- Temporary gingival margin recession resulting from the mechanical compression applied during the procedure;
- Localized pain of short duration after the procedure.
- In rare cases, hypersensitivity reactions related to the materials composing the product may occur. If symptoms persist or worsen, discontinue use of the product and adopt the appropriate clinical management.

WARNING:

- Improper use of POTENZA RETRACTTA may cause gingival trauma, excessive bleeding, pain, local inflammation, or temporary gingival margin recession.


- Do not apply excessive force during insertion of the cord into the gingival sulcus. The application must be performed gently and in a controlled manner, respecting the patient's anatomy and clinical condition.
- Correctly select the cord diameter according to the depth of the gingival sulcus and the patient's periodontal condition. The use of an inappropriate thickness may compromise the effectiveness of gingival retraction or cause tissue injury.
- Do not use the cord on injured, inflamed, infected tissues or in the presence of uncontrolled active bleeding, unless under the clinical evaluation and indication of the dentist.
- When the cord is used in combination with a hemostatic or astringent solution, ensure that the substance is indicated for dental use and strictly follow the instructions for use, contraindications, and warnings provided by the solution manufacturer.
- Do not reuse POTENZA RETRACTTA. The product is intended for single use, and its reprocessing or reuse may compromise its performance and increase the risk of cross-contamination.

INSTRUCTIONS FOR USE:

Step 1: Select the thickness of POTENZA RETRACTTA, identified by the corresponding color, according to the width and depth of the gingival sulcus, as well as the anatomy of the tissue.

Step 2: Whenever necessary, choose a diameter that allows proper accommodation within the gingival sulcus, providing effective retraction without applying excessive pressure to the tissues.

Step 3: Cut the POTENZA RETRACTTA using sharp scissors to a length sufficient to completely encircle the tooth or the area to be treated.

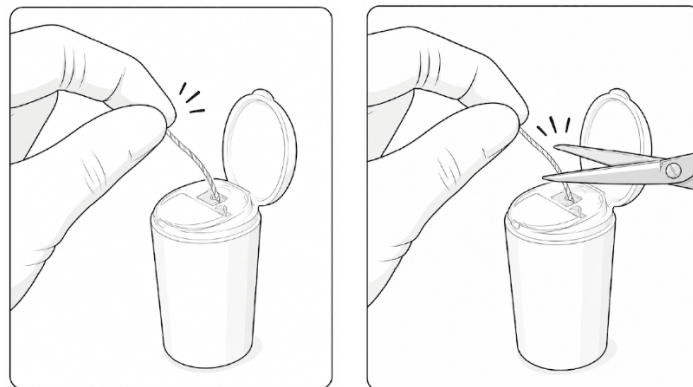Step 4: If necessary, dry the operative field to facilitate visualization and insertion of the retraction cord.

Step 5: POTENZA RETRACTTA may be used dry, providing mechanical retraction of the gingival tissue, or previously soaked in a hemostatic or astringent solution when there is a need to enhance the control of bleeding and sulcular fluids.

Step 6: Using an appropriate retraction cord placement instrument without a cutting edge, carefully insert POTENZA RETRACTTA into the gingival sulcus.

Step 7: Begin the insertion preferably in the interproximal region, where the gingival sulcus tends to be deeper, positioning the cord continuously and uniformly around the entire tooth.

Step 8: Keep the cord in position for the required amount of time. The duration may vary depending on the technique used, gingival anatomy, and whether or not a hemostatic or astringent solution is used.

Step 9: Before removal, if necessary, lightly moisten the cord with water to facilitate its removal and minimize the risk of trauma to the gingival tissue.

Step 10: Carefully remove POTENZA RETRACTTA and proceed with the clinical procedure.

DOUBLE-CORD TECHNIQUE (WHEN APPLICABLE)

For impression procedures, the double-cord technique may be employed:

Step 1: Initially insert a smaller-diameter cord (e.g., #000 or #00) into the bottom of the gingival sulcus to assist in controlling bleeding and sulcular fluids.

Step 2: Next, position a second, larger-diameter cord (e.g., #0 or #1) over the first, promoting lateral displacement of the gingival tissue.

Step 3: Wait approximately 5 minutes, or for the time necessary to obtain adequate gingival retraction.

Step 4: Immediately before taking the impression, remove only the upper cord, keeping the lower cord positioned within the gingival sulcus.

Step 5: Perform the impression according to the technique and impression material selected by the professional.

Step 6: After completing the procedure, carefully remove the remaining cord.

CONSERVATION AND STORAGE:

- Keep POTENZA RETRACTTA in its original, properly sealed packaging until the time of use.
- Store the product in a clean, dry, well-ventilated place, protected from moisture, dust, direct sunlight, and excessive heat sources.
- Store the product at a temperature between 32 °F and 86 °F (0 °C and 30 °C).
- Avoid contact with water, excessive moisture, vapors, or other conditions that may compromise the integrity of the product or its packaging.
- Do not store the product together with chemicals, solvents, or materials that may contaminate or deteriorate its characteristics.
- Do not use the product if the packaging is damaged, punctured, wet, or shows signs of contamination.

TRANSPORT CONDITIONS:

- Transport POTENZA RETRACTTA in its original packaging, preserving its physical integrity and protecting it against impacts, punctures, excessive compression, and moisture.
- Transport the product in a clean, dry place, protected from rain, direct sunlight, and excessive heat sources.
- Transport the product at a temperature between 32 °F and 86 °F (0 °C and 30 °C).
- Avoid transportation conditions that may damage the packaging or compromise the integrity of the product.

DATE OF MANUFACTURE: SEE PACKAGE.

SHELF LIFE: SEE PACKAGE.

PRODUCT FOR EXCLUSIVE USE IN DENTAL OFFICES. FOR USE BY DENTAL PROFESSIONALS ONLY.

MAKE SURE THESE INSTRUCTIONS FOR USE CORRESPOND TO THE PRODUCT TO BE USED.

TO OBTAIN A PRINTED VERSION OF THESE INSTRUCTIONS FOR USE, PLEASE CONTACT THE COMPANY'S TECHNICAL SUPPORT.

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